English	3
Español	13
Français	23

EN Index

Safety Instructions.....	4
Meet no!no!	5
Know Before You no!no!	6
Hot Blades.....	6
Putting no!no! Together.....	7
Working With no!no!.....	7
Hot Blade Care.....	9
Troubleshooting.....	10
Technical Information.....	11

Safety Instructions

READ ALL INSTRUCTIONS BEFORE USING.
ALL WARNINGS AND CAUTIONS SHOULD BE STRICTLY FOLLOWED

As with any electrical devices, certain precautions must be taken in order to ensure your safety.

WARNING – To reduce the risk of burns, fire, electric shock or injury:

- Do not reach for a device that has fallen into water.
- Do not place or store device near a sink or bathtub.
- Do not place in or near water or other liquids.
- Always unplug the power cord from the electrical outlet immediately after use and before cleaning or maintenance.
- Only use this device for its intended purpose as described in this manual.
- This device should never be left unattended when plugged in.
- Close supervision is necessary when this appliance is operated by users under 18 years of age.
- Never operate this device if the cord or plug is damaged in any way.
- This device should only be used with the supplied transformer and cord.
- Keep out of reach of children.
- The thermodynamic wire gets red hot during use; never allow it to touch your skin.
- Do not use no!no! near or rest it on flammable objects such as paper or cloth, etc.

CAUTION

- Do not use this appliance on sunburned skin, open cuts, bruises, blisters, moles or on raised areas such as scars.
- If you have any skin condition other than those listed above, please consult your physician before using no!no!
- Always turn no!no! off with the power switch and then unplug it when you are finished working with it.
- Do not use this appliance on the face, breasts or genitals.
- If you are pregnant, please consult your physician before using no!no!
- In case of burns, apply first aid treatment.
- Before using no!no!, make sure skin surface is dry: No deodorants, No lubricants, No liquids, No creams, No gels, No water - should be applied to skin surface.

Meet **no!no!**™

Congratulations, you are about to enter the world of painless, professional hair removal. As you do, it's important to remember that no!no! is a professional treatment for hair density reduction. These results take time, patience and a strong commitment to the treatment's unique procedures.

Experience shows that hair density reduction results may appear within a few months of beginning treatment. Please keep in mind that each person has their own unique pattern of hair growth and biology and each will experience results at their own pace. These results take time, patience and a strong commitment to the treatment's unique procedures.

Since no!no! is a completely new concept in hair removal, and one you're probably not familiar with, it is important to thoroughly read the user manual and watch the movie before using it for the first time. Be very careful to master the correct use of the device before using it on rounded or sensitive areas.


Do Not use no!no! on the face, breasts or genitals. Do not begin to use no!no! before carefully **READING AND WATCHING** all operating instructions.


Buffer

Hot Blade for long hair

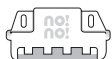
Hot Blade for stubble

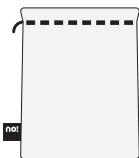
Storage Bag

Cleaning brush

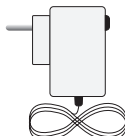Power supply

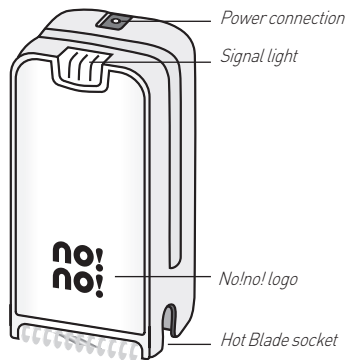Power connection

Signal light



No!no! logo

Hot Blade socket

No!no! device



Know Before You no!no!™

- Some people experience a mild tingling sensation when first starting to no!no!, especially around sensitive areas, this should disappear after the first few treatments.
- Carefully examine the treatment area before the first time you no!no! and before each subsequent treatment. As you track your progress, you should notice that the hair in the treatment area is becoming less dense. These initial results may begin to appear within the first several months.
- For best results use no!no! 2 -3 times a week. no!no! can be used on stubble, it is not necessary to wait for hair to grow long.
- Master no!no! on flatter, less sensitive areas like the legs or forearms. Wait until you have more experience before moving on to more sensitive areas like the bikini line or rounded areas like elbows, knees and underarms.
- Do not use no!no! on sunburned skin, open cuts, bruises, blisters, moles or on raised areas such as scars.
- Those with sensitive skin may experience a burning sensation or find that their skin becomes uncomfortably hot during treatment. To avoid this, take short breaks after repeated gliding on any one area.
- Always make sure thermodynamic wire is completely intact before use. If it is bent or torn, replace the Hot Blade before using it. (fig 4) (fig 5)
- To avoid possible malfunction or harm, make sure to gently push the Hot Blade all the way into the socket. (fig 10).
- Check and clean no!no! frequently especially when working on long hair.

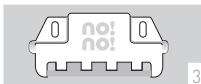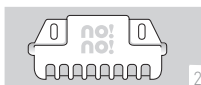WARNING!

Always check the Hot Blade before use. Never use a Hot Blade if the thermodynamic wire is bent, torn or dirty. (fig 4 & 5)

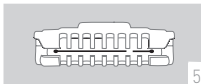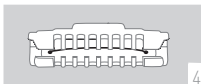Hot Blades

no!no! comes with 2 types of Hot Blades, one for long hair (fig 2) and another for stubble (fig 3). Each Hot Blade has a thermodynamic wire for removing hair. The main difference between them is that the Hot Blade for stubble has a metal strip for more efficient stubble removal.

For optimal treatment results, it is important to use the correct Hot Blade for the length of hair you are working on. The Hot Blade for long hair should be used in the initial stages and whenever hair is approximately 1mm or more. The Hot Blade for stubble should be used on any hair that measures approximately 1mm or less.

Putting **no!no!**™ Together





Remove the safety seal from the no!no! Hot Blade socket (fig 6).

Pick up the device and turn it over. Looking inside the Hot Blade Socket you will see an LED with a metal post on each side. (fig 7) To install the Hot Blade grasp it by both ends.

(fig 8) Line up the two small holes on the Hot Blade with the posts inside the Hot Blade socket and push it gently but firmly into place. Do not use force. To remove the Hot Blade, grasp it by both ends and gently pull it out of the Hot Blade socket.

To install the Hot Blade, hold it so that the logo on the Hot Blade is facing the same direction as the logo on the no!no! (fig 9). Make sure the Hot Blade is seated all the way into the Hot Blade socket. (fig 10)

Next connect the power supply to no!no! (fig 11)

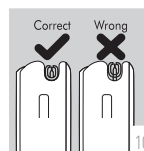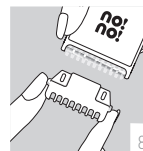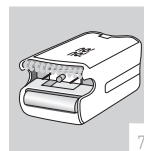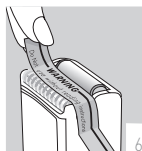Working With **no!no!**™

Skin Preparation

no!no! must be used on clean, dry skin. Always be sure to clean your skin of any gels, oils, lotions or deodorants before you begin to no!no!

Watch the Light


The signal light (fig 12) will guide you throughout your no!no! hair removal experience. It lets you know when you are working at the right speed and warns you if the Hot Blade is torn or missing. The signal light will not be active unless no!no! is in motion.

While gliding no!no! on the skin, the Signal Light should always be facing you. When working at the right speed, the light will be a steady green. If the light flickers or turns off, try adjusting your speed until the light remains a steady green as you glide.

If, at any time, the Signal Light turns red remove the Hot Blade and check it, if the wire is damaged – replace the Hot Blade. Never use a damaged Hot Blade.

Working With no!no!™

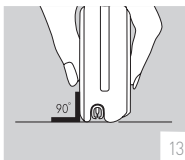
Mastering the Slow Glide

The slow glide technique is crucial to the treatment. To master it, hold no!no! at a 90° angle with the rollers flat against the skin. (pic 13) Then gently glide it over each section you want to treat using a slow, smooth, even motion always pulling it towards your body. Continue to glide over the section you are treating until you feel it is smooth or nearly smooth. Remember: Always glide upwards toward your body in one direction only. And always at a gentle even pace.

Mastering the Buffer

To remove all signs of crystallization and to maximize smoothness, after you complete an area and before moving on to the next section, always use the Buffer. Buffering exfoliates the skin, removes treated hair and gives you an extra smooth surface. To Buff, hold the Buffer flaps (fig 14) and using either a circular or back and forth motion, apply a strong rubbing action to the skin. Always keep Buffers on a clean, dry surface within easy reach whenever you no!no!

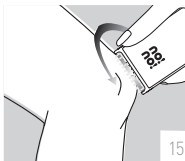For Advanced Users

Because of curves and angles, rounded areas like elbows, knees, wrists, underarms and the bikini line require special care and experience.

Once you have mastered the no!no! technique on legs or forearms, that is you have used no!no! several times and are able to keep the signal light a steady green while gliding and are seeing good hair removal results, you may be ready to try no!no! on more difficult areas.

To work on rounded areas, (figs15) go slowly and allow no!no! to follow the natural curve of the body. (for more information, see the training CD for examples on no!no!ing on rounded areas)

WARNING! Never use no!no! on the face, breasts or genitals.

Be aware that the bikini line can be more temperature sensitive than other areas. If the treatment area becomes uncomfortably warm during treatment, move away from it and work on a different area for awhile. Do not go further than the bikini line.

Treatment

To begin treatment, push the power button on the power supply (fig 16)
When doing the no!no! treatment there are just a few simple rules to remember:

1. Install the appropriate Hot Blade for the length of hair you are working on.
2. Work on one small section of the treatment area at a time, completing one before moving on to the next. Go over each area a few times until you feel that the surface is smooth or nearly smooth.
3. Glide in one direction only – always pulling no!no! upward towards your body with the signal light facing you. (fig 17)
4. Between each section, use the Buffer.
5. Periodically check the Hot Blade for hair debris during treatment and clean it whenever necessary.
6. After you finish the entire treatment area, if you started with the Hot Blade for long hair, switch to the Hot Blade for stubble; the one with the metal strip, and repeat.

Hot Blade Care


Hot Blades are designed to last 2-3 hours of active use. This includes only the time when the Hot Blade is in actual working contact with the skin. The actual life of a Hot Blade, however, depends on several varying factors among them; frequency of use and the area and type of hair being treated.

It is important to keep Hot Blades free of hair debris. Check and clean them periodically during each treatment session. Allowing hair to build up in the Hot Blade may cause it to burn out prematurely or malfunction. To clean the Hot Blade, unplug or turn off your no!no!, remove the Hot Blade from its socket and gently brush it clean with the cleaning brush supplied with your no!no! (pic 18)

Replace the Hot Blade:

- if you begin to notice that it is not removing hair as effectively as before
- if the wire is damaged in any way
- if the signal light turns red

For instructions on replacing Hot Blades, see page 7.

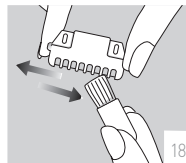Troubleshooting

If no!no! does not start when you glide it on the skin:

Make sure the On/Off button is switched to On

Make sure that the power cord is plugged firmly into the device and the electrical wall outlet.

Make sure the Hot Blade is intact. Replace the Hot Blade if necessary.

If the thermodynamic wire is severed (pic) replace the Hot Blade

If Thermodynamic Wire is bent (pic) replace the Hot Blade

If the Signal light indicator turns red while you glide no!no! on your skin:

Make sure Hot Blade is correctly installed in Hot Blade socket

Make sure Hot Blade is not torn.

If the device is not working properly, or if it has been damaged or dropped into water, return it to a designated service center for examination or repair.

Signal lights	Green	Red
Correct gliding speed	On and steady	Off
No Hot Blade in socket	Off	On
Damaged Hot Blade	Off	On
Incorrect gliding speed	Flickering	Off
No movement	Off	Off

Technical Information

TECHNICAL SPECIFICATIONS:

Weight: 130gr

Physical Dimensions: 38 X 33 X 81mm
(1.5 X 1.3 X 3.2 in)

Mode of Operation: Continuous

OPERATING CONDITIONS:

Temperature: 5° C - 30° C

Humidity: 30% - 80% RH

Rating: 100-240VAC
50-60 HZ
0.3W

TRANSPORT AND STORAGE

Temperature: -20° C to 80° C

Humidity: 05-95% RH
0.1-1.5 atm.

Atmospheric pressure: 100-240VAC
50-60 HZ, 0.3W

ELECTRICAL SAFETY:

Class II equipment - double insulation

 power supply only)

This appliance is not suitable for use in the presence of flammable anesthetic mixture with air or with Oxygen or Nitrous Oxide.

The symbol  on the label affixed to this appliance means "Attention, consult accompanying documents".

The Quality Management System complies with the Quality Management Standard ISO 13485.

no/No! complies with IEC 60335 and with Low Voltage Directive (Electric Safety) LVD: 2006/95/EC

The design and technical specifications are subject to change without notice.

Manufacturer:

40 Ramland Road South, Suite 10
Orangeburg, New York 10962 USA

EC Representative:

Obelis S.A Av.de Tervuren 34 Bte. 44
B-1040 Brussels, BELGIUM

Tel: 32 (0) 2 732 59 54

Fax: 32 (0) 2 732 60 03

GSM: 07545 46 60

Email: obelis@info.be

Service Centers:

For service, please refer to the warranty to find the nearest service center or go to our website

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